



F-Scan™ GO

Upgrading from F-Scan Versatek™ systems to F-Scan™ GO

What is F-Scan GO?

In-Shoe Foot function & Gait Analysis System

F-Scan GO is the next generation of the F-Scan In-Shoe Pressure Measurement system. The F-Scan In-Shoe system provides dynamic pressure, force and timing information for foot function and gait analysis. Information obtained from the F-Scan is used in real-world applications, like designing and testing orthotics, offloading diabetic feet and evaluating footwear and techniques in elite athletes.

The new F-Scan GO system features completely cordless electronics and sampling rates up to 500 Hz over Wi-Fi, with additional onboard SD card storage for extended range.



New F-Scan GO



Legacy F-Scan Systems




Why Upgrade to F-Scan GO?

- 1 **Complete freedom of movement that preserves natural gait**
- 2 **No cables or cords**
- 3 **Fast setup** – connect sensors and “GO”!
- 4 **Extended recording range (datalogging) and 32GB onboard storage**
- 5 **FootView software** – improved interface, more intuitive, and enhanced for use with a tablet.


Best-In-Class Sensors

Tekscan sensors are the thinnest, and the highest-resolution sensor on the market today. Only Tekscan offers trimmable sensors which can be perfectly sized for an unhindered fit of any shoe or foot, ensuring that data capture is complete and accurate.

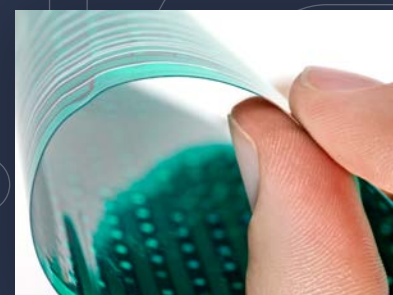


FootVIEW Software

Modernized interface streamlines data collection and analysis, especially for the Biomechanics research user.

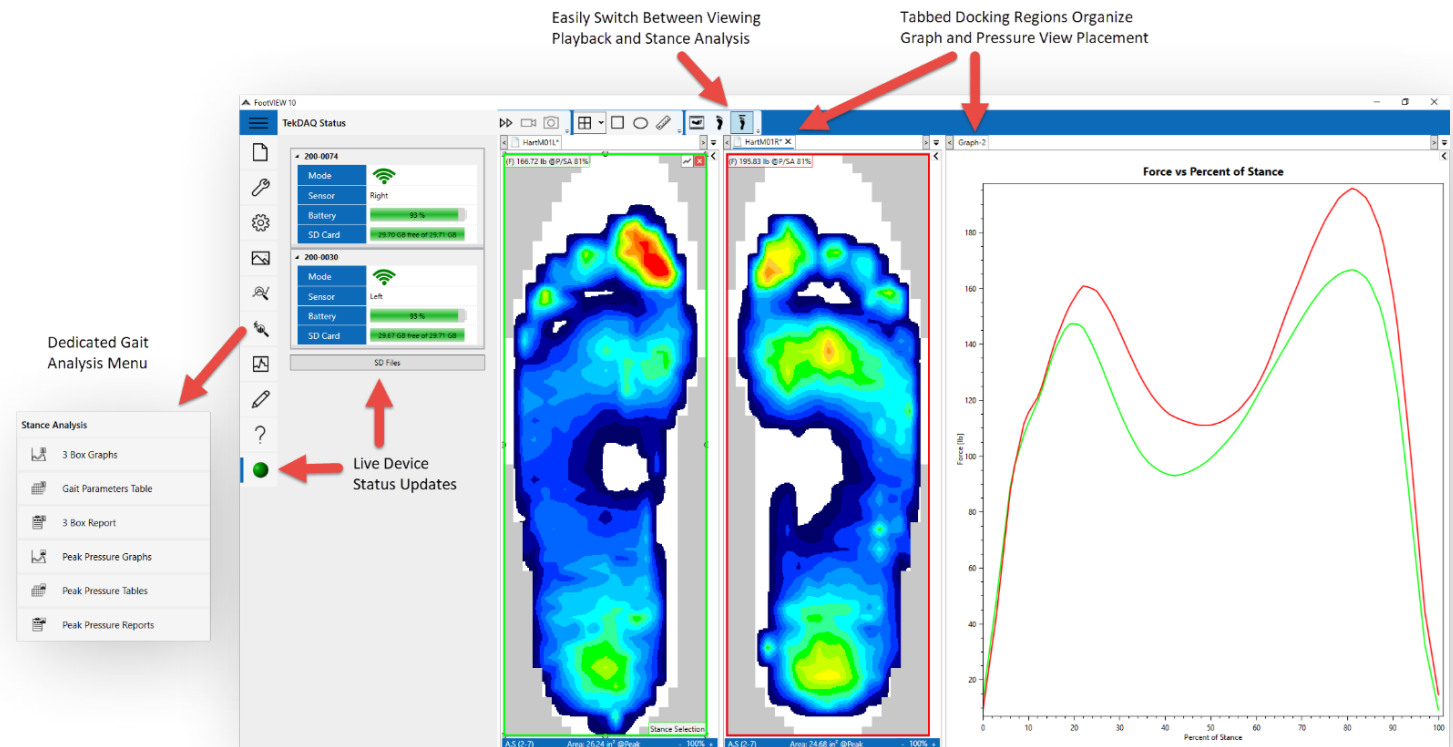


Peak Pressure Analysis

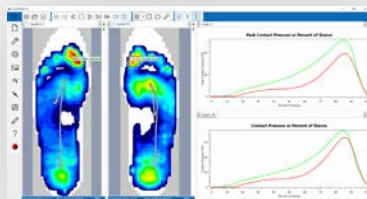


- Identifies and quantifies the peak pressure areas
- Confirms the efficacy of offloading treatments
- Generates a report showing before and after comparison

3-Box Analysis

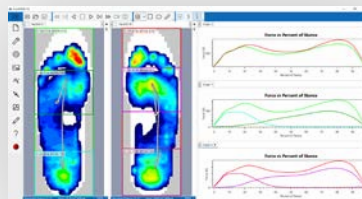

- Automatically segments the foot into heel and forefoot regions
- Easily display force curves to measure symmetry and timing of foot function
- Export to customizable, printable report templates for documentation and patient education

Automated Analysis with Reports

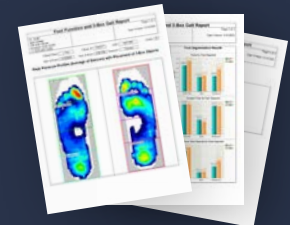

- Center of Force Trajectory
- Real-time status view of TekDAQ devices
- Wi-Fi enabled communication and control for TekDAQ devices
- ASCII Export of pressure and graph data for analysis
- Optional Data Reader Toolkit library for direct integration with external software for analysis of recorded data
- Synchronization with external devices through external trigger input controls

Key Applications



Research

Biomechanics researchers have come to trust Tekscan to deliver accurate and reliable in-shoe force and pressure measurement systems and data to suit a limitless range of dynamic applications. Countless research-validated studies support this fact; and the development of F-Scan GO continues Tekscan's position as a true leader in the field of in-shoe gait analysis.



Gait Analysis and Lower Extremity Function

F-Scan GO is a useful tool for diagnosing pathologies; and evaluating the effect of offloading treatments, footwear modifications, and physiological changes, such as joint manipulations. In many cases, it is an invaluable aid in treating patients as they search for solutions to help alleviate their pain.




Orthopedics, Physical Therapy and Rehabilitation

An objective, reliable tool for PT and Rehab professionals to measure patient progress and help them understand their pathologies and treatment. Additionally, F-Scan GO can differentiate a practice and increase revenue.



Wound Care

Tekscan pressure analysis systems are used extensively by pedorthists, podiatrists and others working in the field of diabetic care. F-Scan GO is a critical tool that assists in the management of in shoe pressure to treat and prevent diabetic ulcers – primarily through the use of designing, manufacturing and modifying orthotics.

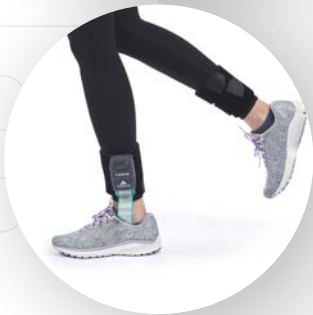

Sports Medicine and Performance

In-shoe gait analysis systems provide an accurate, objective way to capture baseline information, evaluate the effectiveness of strength and conditioning programs, track progress over time and develop injury prevention programs. F-Scan GO is ideal for these applications as it offers complete freedom of movement, a faster scan rate (500 Hz), and the ability for data logging.

New In-Shoe Sensors

New Model

3010 F-Scan In-Shoe Sensors

Slightly thicker substrate to aid sensor durability, while maintaining “thinnest sensor on the market”

E-Tab connector – thinner profile for a more streamlined connection

New routing layout reduces chances of disconnections

Technical Specifications



The F-Scan GO system consists of:

- (2) F-Scan GO TekDAQ wearable Wi-Fi data acquisition units
- (2) Removable rechargeable battery packs
- (2) Battery Chargers
- (2) 32GB formatted Micro-SD cards
- (2) Velcro ankle bands
- (20) Trimmable model 3010 F-Scan GO sensors (10 pairs)

Sensor		Electronics	
Max. # of Sensels	Up to 966	Included	2 TekDAQ Wi-Fi Data Collection Units
Sensor Thickness	0.28 mm	Max Scan Rate	500 Hz
Resolution	3.9 per cm2 / 25 sensels™ per in2	Weight	266g / 9.4 oz (133g / 4.7 oz per leg)
Pressure Rating	125 PSI / 862 kPa	Onboard Storage	32 GB – Expandable with Micro SD
Sensitivity Adjustment	+/- 3x Adjustment of Pressure Rating	Battery Life	2 Hours
Sensor Sizes	Trimmable Sensors: Ranging from Men's 14 (USA) - Smallest size Girls Size 2	Software	FootVIEW



617.464.4282 1.800.248.3669 info@tekscan.com

Ask us about special trade-in offers for your existing F-Scan!